

Equipamientos Analíticos S.A

DESDE 1989 | BRINDANDO SOLUCIONES

FAMILY THREE

High-Pressure Pumps

The Floxlab signature line — high-pressure syringe pumps, pressure controllers and automatic valves engineered for the most demanding fluid-handling applications in research and industry.

HIGH-PRESSURE PUMPS

Syringe Pumps · Continuous Flow · Pressure Controllers · H₂ / CO₂ Solutions

40,000
psi
Max pressure

160
°C
Heating option

Pulse-free
Flow output

SS / Hast.
Wetted parts



OVERVIEW

Pump Product Range

A complete ecosystem for every flow and pressure challenge

01

SYRINGE PUMPS

High-pressure precision flow

CF · BTSP / BTDP · BFSP / BFDP · PV

02

PRESSURE CONTROLLERS

Confining & back-pressure

ACP · ACP-150 · ABPR

SYRINGE PUMPS

CF-series

Continuous Flow Pump — pulse-free, unlimited flow output



CF pump (front view)



CF pump (rear view)

DUAL-PISTON WORKING PRINCIPLE

The Floxlab CF Series delivers unparalleled accuracy and repeatability, providing critical pulse-free, high-pressure flow metering. Dual motor-driven pistons employ a unique pre-pressurization step — one cylinder dispenses fluid while its counterpart rapidly refills and pre-pressurizes, ensuring a seamless, pulse-free transition.

CF-SERIES SPECIFICATIONS

Model	Pressure (psi)	Volume (ml)	Max Flow (ml/min)
CF 3	3,000	2 × 40	80
CF 6	6,000	2 × 32	55
CF 12	12,000	2 × 15	30
CF 15	15,000	2 × 12	25
CF 20	20,000	2 × 10	15
CF 30	30,000	2 × 5	5

APPLICATIONS

Core flooding · EOR · Chromatography · Analytical dispensing · Process calibration

SYRINGE PUMPS

BTSP-series

Benchtop High-Pressure Single Syringe Pumps



BTSP pump



BTSP with cooling jacket



BTSP with heating mantle

The BTSP series provides ultra-precise, pulse-free fluid flow control for both inorganic fluids (e.g., brine, solvents) and organic fluids (e.g., hydrocarbons, microbial solutions). It delivers controlled pressure, flow rate, and volume with capacities up to 40,000 psi, and supports both constant-pressure and constant-flow modes. A high-precision pressure transducer enables real-time monitoring, while two manual isolation valves regulate inlet and outlet flow. The system includes data acquisition and control software for monitoring and analysis.

STANDARD MODELS

Model	Vol (ml)	P (psi)	Q (ml/min)
BTSP 20-40	20	40,000	7.5
BTSP 50-30	50	30,000	20
BTSP 125-20	125	20,000	30
BTSP 175-15	175	15,000	30
BTSP 250-10	250	10,000	60
BTSP 500-5	500	5,000	130
BTSP 1000	1,000	1,875	250

EXTENDED MODELS

Model	Vol (ml)	P (psi)	Q (ml/min)
BTSP 150-30	150	30,000	20
BTSP 250-20	250	20,000	35
BTSP 300-15	300	15,000	40
BTSP 500-10	500	10,000	70
BTSP 1000-5	1,000	5,000	130

OPTIONS

Touch-screen · Heating mantle · Cooling jacket · Mixer · Hastelloy

SYRINGE PUMPS

BTDP-series

Benchtop High-Pressure Dual Syringe Pumps



BTDP pump



BTSP with Heating mantles



BTDP with mixers

The BTDP series combines two BTSP pump modules, enhancing system capabilities. Equipped with a pair of pneumatic powered three-way valves for both tank inlet and outlet delivery, the dual pump provides uninterrupted, continuous fluid flow across a comprehensive range of both pressure and temperature conditions. All functionalities and optional features inherent to the BTSP series pumps remain accessible as both individual pump modules can still operate autonomously.

STANDARD MODELS

Model	Vol (ml)	P (psi)	Q (ml/min)
BTDP 20-40	2 x 20	40,000	7.5
BTDP 50-30	2 x 50	30,000	20
BTDP 125-20	2 x 125	20,000	30
BTDP 175-15	2 x 175	15,000	30
BTDP 250-10	2 x 250	10,000	60
BTDP 500-5	2 x 500	5,000	130
BTDP 1000	2 x 1,000	1,875	250

EXTENDED MODELS

Model	Vol (ml)	P (psi)	Q (ml/min)
BTDP 150-30	2 x 150	30,000	20
BTDP 250-20	2 x 250	20,000	35
BTDP 300-15	2 x 300	15,000	40
BTDP 500-10	2 x 500	10,000	70
BTDP 1000-5	2 x 1,000	5,000	130

BTDP = 2 x BTSP with 3-way valves **OPTIONS**

Touch-screen · Heating mantle · Cooling jacket · Mixer · Hastelloy



ENERGY TRANSITION

H₂ & CO₂ Liquid Pumps

Specialized syringe pumps for hydrogen and liquid CO₂ applications

FloXlab BTSP pumps can be specially configured for hydrogen (H₂) and liquid CO₂ service — critical fluids for green hydrogen research, CCS/CCUS (carbon capture & storage), and enhanced oil recovery. Custom thermal jackets, Hastelloy wetted parts, and high-pressure ratings guarantee safe and accurate handling of these challenging fluids.

H₂ — HYDROGEN SERVICE



BTSP 500-5 H₂ — Hydrogen application

- Pressure up to 10,000+ psi
- Specialized seals & gaskets
- SS 316 wetted materials
- Low-permeation design
- Safety interlocks & leak detection

Green H₂ · Fuel cells · Electrolyzers · Storage

CO₂ — LIQUID CO₂ SERVICE



BTSP CO₂ pump with cooling bath

- Cooling jacket (sub-ambient)
- Liquid CO₂ state maintained
- Supercritical CO₂ (scCO₂) capable
- Thermally insulated cylinder
- High-pressure (up to 10,000 psi)

CCS / CCUS · scCO₂ chemistry · CO₂ EOR · Storage

SYRINGE PUMPS

BFSP / BFDP-series

Bench-Floor Syringe Pumps — large volume & mobility



BFSP / BFDP pumps

BENCH-FLOOR SYRINGE PUMPS

LARGE VOLUME

Bench-floor configurations for applications requiring large volumetric capacities (500 to 1,000 ml) with pressure up to 25,000 psi.

KEY FEATURES

- Sturdy chassis with 4 heavy-duty casters
- BFDP dual version with 3-way valves
- Touch-screen control panel
- Heating mantle up to 150 °C
- SS or Hastelloy wetted parts

BFSP / BFDP SPECIFICATIONS

Model	Volume (ml)	Pressure (psi)	Max Flow (ml/min)
BFSP 500-15	500	15,000	40
BFSP 1000-15	1,000	15,000	80
BFSP 500-25	500	25,000	50
BFDP 500-15	2 × 500	15,000	40
BFDP 1000-15	2 × 1,000	15,000	80
BFDP 500-25	2 × 500	25,000	50

COMMON FEATURES

Shared Features

Common across all Floxlab syringe and continuous-flow pumps



Touch-screen control panel — computer-independent operation



Fluid stirrer

01

FLEXIBLE CONTROL

Modes: time, pressure, flow rate, volume, or trigger events

02

TOUCH-SCREEN PANEL

Intuitive front-panel interface for computer-independent operation

03

DELIVER OR RECEIVE

Bidirectional operation: deliver or receive fluid with equal precision

04

FLUID STIRRER

Optional integrated stirrer for homogenized fluid handling

05

WETTED MATERIALS

High-grade stainless steel standard · Hastelloy for corrosive fluids

06

THERMAL OPTIONS

Heating mantle up to 150 °C · Cooling jacket for sub-ambient Temp.

EVERY FLOXLAB PUMP: Ready-to-use · Data acquisition software · System supervision · Advanced analytics


SPECIALTY PUMPS

PV-series & Pharmacia

Pressure Volume pumps and liquid metering solutions

PV-SERIES

Pressure Volume Pump

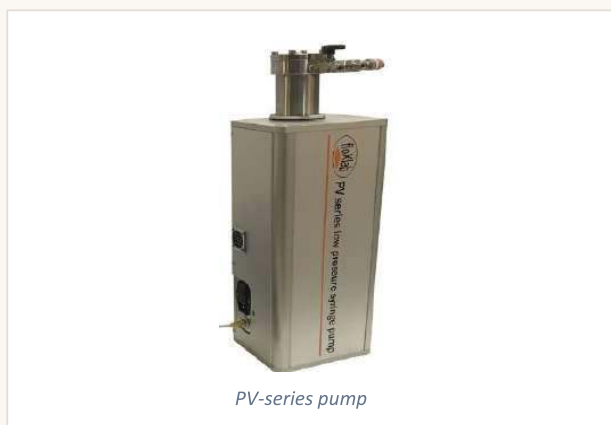

FEATURES


- Single cylinder configuration
- Ambient temperature operation
- Stainless steel wetted parts
- Supervision software included
- Ethernet connection for data logging
- Horizontal / vertical position

Model	P (psi)	Vol (ml)	Q (ml/min)
PV 200	600	200	70
PV 500	225	500	150

Water pore volume change · Pressure control · Soil studies

PHARMACIA

Liquid Metering Pump

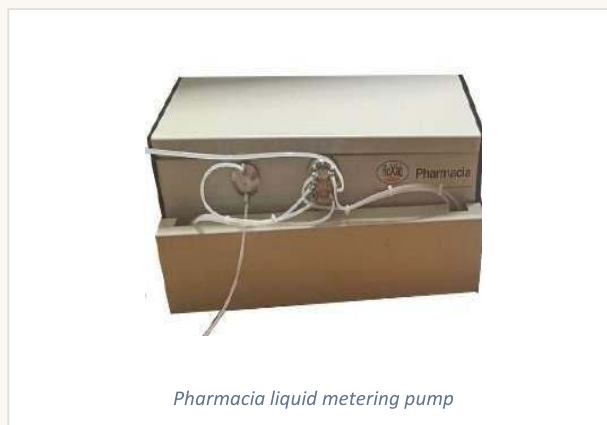

FEATURES

- Dual-piston antiphase architecture
- Quasi-continuous flow output
- Automatic switchover valve
- Constant flow OR constant pressure modes
- Stainless steel wetted parts
- Computer control via dedicated software

Model	P (psi)	Vol (ml)	Q (ml/min)
Pharmacia	300	2 × 15	20

Solvent injection · Liquid continuous flow



PRESSURE CONTROLLERS

ACP · ACP-150 · ABPR

Dedicated systems for confining, pore and back-pressure regulation

ACP-series

Automated confining pressure

- Servo-controlled confining
- Compatible with triaxial cells
- Closed-loop PID control
- Pressure ramp & hold modes
- High accuracy & repeatability



ACP-150 series

Pressure & temperature controller

- Combined P + T regulation
- Operating temp up to 150 °C
- Reservoir condition simulation
- Heating mantle with PID
- Safety interlocks included



ABPR-series

Automated back-pressure regulator

- Precise pore/outlet pressure
- Maintains downstream conditions
- Pneumatic dome BPR design
- Pressure 10,000 to 20,000 psi
- Gas and liquid service



APPLICATIONS

FloXlab pumps power a wide range of laboratory and industrial systems



H₂ PERMEAMETER
Hydrogen flow studies



CHEMISTRY & PROCESS
Heavy feedstock Pilot plant



FRACLAB
Hydraulic fracturing

SECTORS: Rock Mechanics · Reservoir Engineering · H₂ & CCUS · Hydraulic Fracturing · Petrophysics · Chemistry & Process

R r

fl fl fl fl fl fl fl

fffi @ P6 R r P (fi xfi fi fiV fiwfi fiPP6 fi P x (V 8 x
fi fi P (xfi xfi P x (fi fi P x V x Xfi (fix fi(Pfi
fiP fi ff @ P6 (fi P x fi fi V fi wfffi fi P6 fi 8 P Pfixw
P(fi fi (V PV (fi fiVfP V wfi fi P x Vi P VP



R- R

3 (ff fi (ff
fi fiV V&fix6PV V



2

s 8 ff(ff fi fi
P6 P x fi



R r

P fi fi(P ff @ P6
fi 8 Pfi

R
· r 2

(P x P ff @ P6
p V fi

R @ P6 (fi fi fi u fi P P fi

fl fl fl fl

fi x w P fi fi @P Vi fi fiH V PVS 8 P6 V 8 P6 V



TECHNICAL BROCHURE

AUTO-VALVE

S E R I E S

Solenoid-actuated pneumatic valves · High-temperature, high-pressure

Up to 30,000 psi

Up to 160 °C

SS / Hastelloy

OVERVIEW

Auto-Valve Series — Key facts

DESIGN HIGHLIGHTS

- Solenoid-actuated pneumatic valves
- Computer-controlled via solenoids
- Corrosion-resistant construction
- No volume or pressure change during cycle
- Two configurations: ON/OFF & 3-way 4-pos

TECHNICAL SPECS

PRESSURE

5,000 – 30,000 psi

TEMPERATURE

80 °C / HT 160 °C

WETTED PARTS

Stainless Steel / Hastelloy

SEALS

Viton (std) · PU on 30K

AIR SUPPLY

80–115 psi pneumatic

KEY BENEFITS

Why engineers choose Floxlab Auto-Valves

Computer-Controlled

Remote, repeatable solenoid actuation.

Pneumatic Actuation

Fast response with 80–115 psi supply.

High Pressure

Proven from 5K up to 30K psi.

High Temperature

HT option rated to 160 °C.

Corrosion Resistant

SS or Hastelloy wetted parts.

Constant Volume

No volume/pressure change in cycle.

2 - WAY VALVES

ON/OFF Valves

Single valve configuration — pressure families from 5K to 30K psi



CONFIGURATION

TYPE

Single valve, ON/OFF

ACTUATION

Pneumatic, solenoid-piloted

AIR SUPPLY

80–115 psi

SEAL

Viton (std) · PU on 30K

APPLICATION

Ideal for isolation, shut-off and on/off control in high-pressure lines. Computer-controlled via solenoids with fast, clean pneumatic actuation.

PRESSURE FAMILIES

5,000 psi

1/4" LOK
80/160 °C
SS / Hast.

10,000 psi

1/8" LP
80/160 °C
SS / Hast.

10K psi HF

1/4" SLP
80/160 °C
High flow

20,000 psi

1/4" HP
80/160 °C
SS / Hast.

30,000 psi

1/4" HP
40 °C · PU
SS / Hast.

COMPLETE MODEL LIST

Model	Pressure	Temp	Wetted	Cv	Int. Vol.	Connection	Dimensions
SS05-AT	5,000 psi	80 °C	SS	0.10	0.50 cc	1/4" LOK	5.1 × 10.6 cm
SS05-HT	5,000 psi	160 °C	SS	0.10	0.50 cc	1/4" LOK	5.1 × 10.6 cm
SH05-AT	5,000 psi	80 °C	Hastelloy	0.10	0.50 cc	1/4" LOK	5.1 × 10.6 cm
SH05-HT	5,000 psi	160 °C	Hastelloy	0.10	0.50 cc	1/4" LOK	5.1 × 10.6 cm
SS10-AT	10,000 psi	80 °C	SS	0.055	0.11 cc	1/8" LP	5.1 × 6.8 cm
SS10-HT	10,000 psi	160 °C	SS	0.055	0.11 cc	1/8" LP	5.1 × 6.8 cm
SH10-AT	10,000 psi	80 °C	Hastelloy	0.055	0.11 cc	1/8" LP	5.1 × 6.8 cm
SH10-HT	10,000 psi	160 °C	Hastelloy	0.055	0.11 cc	1/8" LP	5.1 × 6.8 cm
SS10HF-AT	10,000 HF	80 °C	SS	0.085	0.44 cc	1/4" SLP	5.7 × 8.9 cm
SS10HF-HT	10,000 HF	160 °C	SS	0.085	0.44 cc	1/4" SLP	5.7 × 8.9 cm
SH10HF-AT	10,000 HF	80 °C	Hastelloy	0.085	0.44 cc	1/4" SLP	5.7 × 8.9 cm
SH10HF-HT	10,000 HF	160 °C	Hastelloy	0.085	0.44 cc	1/4" SLP	5.7 × 8.9 cm
SS20-AT	20,000 psi	80 °C	SS	0.05	0.22 cc	1/4" HP	5.7 × 8.0 cm
SS20-HT	20,000 psi	160 °C	SS	0.05	0.22 cc	1/4" HP	5.7 × 8.0 cm
SH20-AT	20,000 psi	80 °C	Hastelloy	0.05	0.22 cc	1/4" HP	5.7 × 8.0 cm
SH20-HT	20,000 psi	160 °C	Hastelloy	0.05	0.22 cc	1/4" HP	5.7 × 8.0 cm
SS30-AT	30,000 psi	40 °C	SS / PU	0.05	0.30 cc	1/4" HP	6.8 × 10.0 cm
SH30-AT	30,000 psi	40 °C	Hast./PU	0.05	0.30 cc	1/4" HP	6.8 × 10.0 cm

3 - WAY VALVES

4-Position Valves

Double-valve assembly for fluid routing and switching applications



CONFIGURATION

TYPE

Double valve, 3-way 4-pos

ACTUATION

Solenoid-piloted pneumatic

AIR SUPPLY

80–115 psi

SEAL

Viton (std) · PU on 30K

APPLICATION

Dual-valve configuration with three ports and four switching positions. Ideal for fluid routing, directional switching and multi-line control.

PRESSURE FAMILIES

5,000 psi

1/4" LOK
80/160 °C
SS / Hast.

10,000 psi

1/8" LP
80/160 °C
SS / Hast.

10K psi HF

1/4" SLP
80/160 °C
High flow

20,000 psi

1/4" HP
80/160 °C
SS / Hast.

30,000 psi

1/4" HP
40 °C · PU
SS / Hast.

COMPLETE MODEL LIST

Model	Pressure	Temp	Wetted	Cv	Int. Vol.	Connection	Dimensions
DS05-AT	5,000 psi	80 °C	SS	0.10	0.95 cc	1/4" LOK	5.1 × 15.9 cm
DS05-HT	5,000 psi	160 °C	SS	0.10	0.95 cc	1/4" LOK	5.1 × 15.9 cm
DH05-AT	5,000 psi	80 °C	Hastelloy	0.10	0.95 cc	1/4" LOK	5.1 × 15.9 cm
DH05-HT	5,000 psi	160 °C	Hastelloy	0.10	0.95 cc	1/4" LOK	5.1 × 15.9 cm
DS10-AT	10,000 psi	80 °C	SS	0.055	0.31 cc	1/8" LP	5.1 × 12.0 cm
DS10-HT	10,000 psi	160 °C	SS	0.055	0.31 cc	1/8" LP	5.1 × 12.0 cm
DH10-AT	10,000 psi	80 °C	Hastelloy	0.055	0.31 cc	1/8" LP	5.1 × 12.0 cm
DH10-HT	10,000 psi	160 °C	Hastelloy	0.055	0.31 cc	1/8" LP	5.1 × 12.0 cm
DS10HF-AT	10,000 HF	80 °C	SS	0.085	0.74 cc	1/4" SLP	5.7 × 15.2 cm
DS10HF-HT	10,000 HF	160 °C	SS	0.085	0.74 cc	1/4" SLP	5.7 × 15.2 cm
DH10HF-AT	10,000 HF	80 °C	Hastelloy	0.085	0.74 cc	1/4" SLP	5.7 × 15.2 cm
DH10HF-HT	10,000 HF	160 °C	Hastelloy	0.085	0.74 cc	1/4" SLP	5.7 × 15.2 cm
DS20-AT	20,000 psi	80 °C	SS	0.05	0.62 cc	1/4" HP	5.7 × 14.4 cm
DS20-HT	20,000 psi	160 °C	SS	0.05	0.62 cc	1/4" HP	5.7 × 14.4 cm
DH20-AT	20,000 psi	80 °C	Hastelloy	0.05	0.62 cc	1/4" HP	5.7 × 14.4 cm
DH20-HT	20,000 psi	160 °C	Hastelloy	0.05	0.62 cc	1/4" HP	5.7 × 14.4 cm
DS30-AT	30,000 psi	40 °C	SS / PU	0.05	0.70 cc	1/4" HP	6.8 × 17.2 cm
DH30-AT	30,000 psi	40 °C	Hast./PU	0.05	0.70 cc	1/4" HP	6.8 × 17.2 cm

Nomenclature: D = Double (3-way 4-pos) · S = Stainless · H = Hastelloy · 05/10/20/30 = pressure × 1,000 psi · HF = High Flow